

3 configurations providing expansion capabilities, adapted to your activity.



AIA-900

- Panel or single assay
- Up to 5 different tests per tube
- “auto mode”: place the samples, the AIA-PACKs & press START

Dimensions: W 890 mm, D 665 mm, H 642 mm
Weight: 122.5 kg

AIA-900 + 9tray Sorter

- 180 tests on board
- Up to 9 different tests per tube
- Smooth switch from “auto mode” to “sorter mode”

Dimensions: W 1290 mm, D 665 mm, H 642 mm
Weight: 163 kg




AIA-900 + 19tray Sorter

- 380 tests on board
- Up to 19 different tests per tube
- Fully automated with “sorter mode”

Dimensions: W 1520 mm, D 665 mm, H 1247 mm
Weight: 273 kg

The new generation of standalone analysers, flexible and fully scalable.

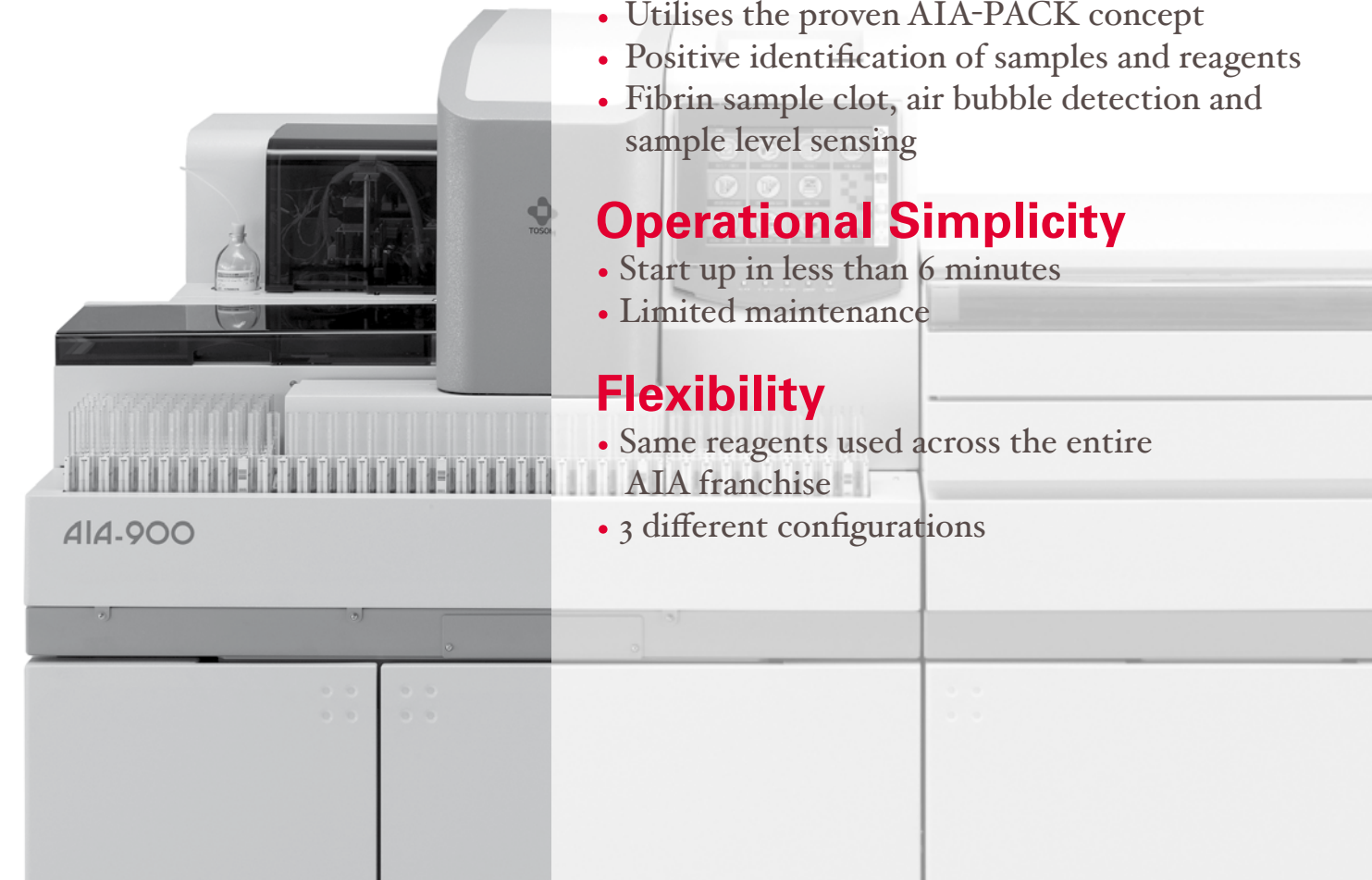

Reliability

- Utilises the proven AIA-PACK concept
- Positive identification of samples and reagents
- Fibrin sample clot, air bubble detection and sample level sensing


Operational Simplicity

- Start up in less than 6 minutes
- Limited maintenance

Flexibility

- Same reagents used across the entire AIA franchise
- 3 different configurations

•AIA-PACK• and •AIA• are the registered trademarks of Tosoh Corporation in EU, etc.

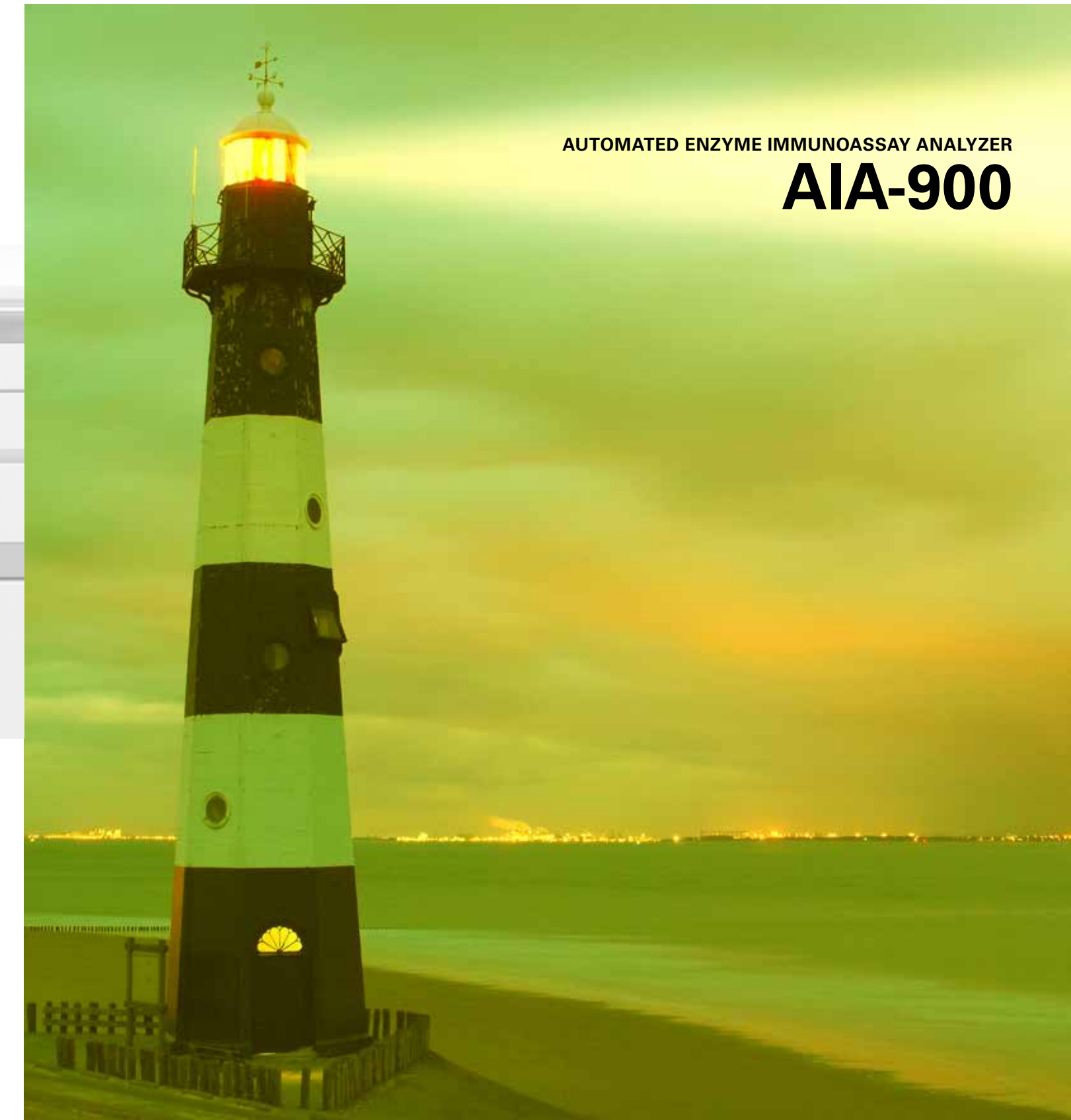

AUTOMATED ENZYME IMMUNOASSAY ANALYZER

AIA-900

TOSOH EUROPE N.V.
Transportstraat 4 - 3980 Tessenderlo - BELGIUM
Tel: +32 (0)13 66 88 30 - Fax: +32 (0)13 66 47 49
www.tosohbioscience.eu

TOSOH BIOSCIENCE LTD
The Business Centre - Edward Street, Redditch
B97 6HA Worcestershire - UNITED KINGDOM
Tel: +44 (0) 1527 592901 - Fax: +44 (0) 1527 592902
www.tosohbioscience.eu

TOSOH CORPORATION - BIOSCIENCE DIVISION
Shiba-Koen First Bldg. 3-8-2, Shiba, Minato-ku,
Tokyo 105-8623 - JAPAN
Tel: +81-3-5427-5181 - Fax: +81-3-5427-5220
www.tosoh.com

revision: AIA900 EN Rev. 28062013.01 - 07/2013 www.hepcasting.be

TOSOH BIOSCIENCE



The AIA-PACK. A unique concept.

- Lyophilised reagents allow long AIA-PACK shelf-life*
- 90 days calibration stability
- Most of the analytes require a 10 minute incubation time, delivering the first result within 20 minutes
- Low risk for contamination or carry-over between samples and reagents:
 - 1 AIA-PACK = 1 test
 - The complete assay (incubation, wash and detection) is carried out in the AIA-PACK
- Disposable pipette tips for samples and reagents
- The single cup per test combined with low volume consumption of reagents reflects Tosoh's environmental awareness

* 12-months from date of manufacture



A ready-to-use "all-in-one" test cup

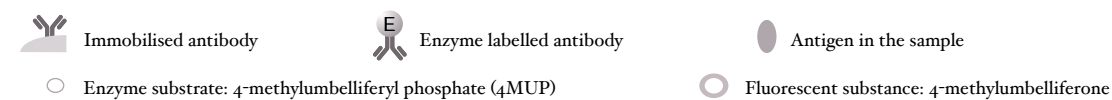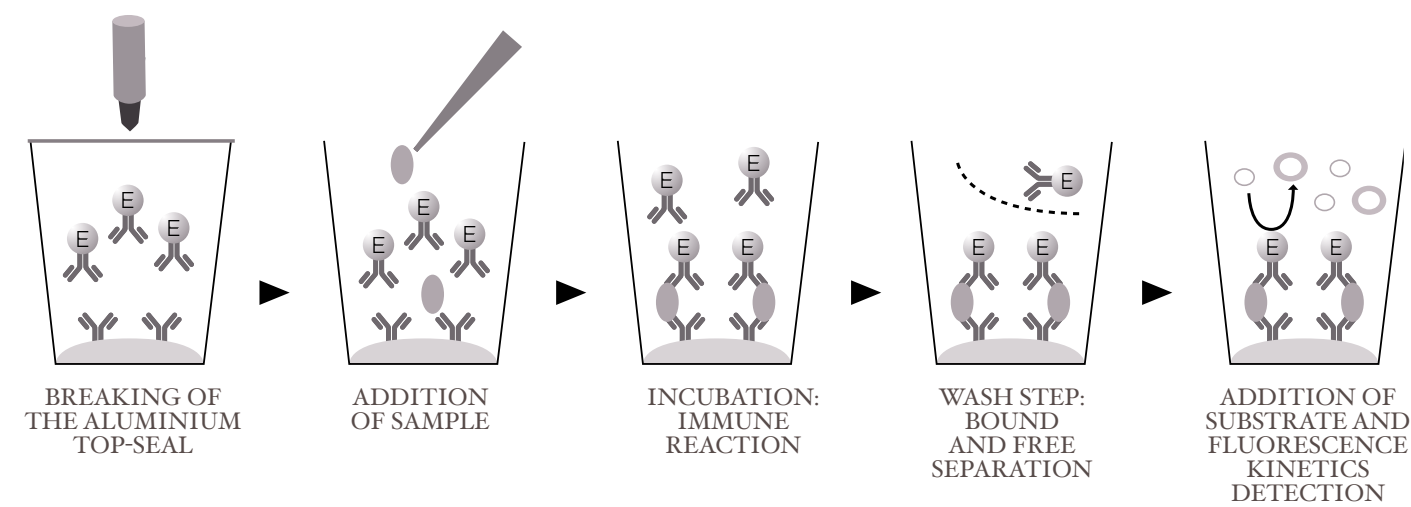

Throughput: up to 90 tests/hour.

Reliability

- Capacitive sample level sensing
- Positive identification of primary tubes
- Absence of carry-over due to the use of disposable pipette tips
- Well established unique AIA-PACK concept



User-friendly operation

- Intuitive touch screen
- Automatic storage of all data on USB memory
- Bidirectional interface to LIS and external QC package
- Limited maintenance
- Start-up in less than 6 minutes



"auto mode"



"sorter mode"

Simple loading

- Load cups and tubes for Loader Model
- Load tubes for Sorter Model

Continuous operation Walk away capability

- Continuous loading of samples, reagents & consumables during operation
- Automated pre-treatment of samples
- Sample loading capacity: 100 (for 9tray & 19tray sorters)
- 2 calibration curves simultaneously active per analyte
- STAT sample rack recognition